

# Oral Syringes and Accessories



## Safe and Accurate Dispensing of Oral Medications

We offer a full-line of oral syringes, bottle adapters, tips, and stands. Our unique oral syringe is a two-piece dispenser that uses a patented wiper design instead of a rubber grommet to form a complete seal between the plunger and the barrel. With no small siliconized rubber parts on the plunger tip, the risk of choking and product contamination is greatly reduced. Custom colored components, color printing, and wrapping options available.



**73 YEARS OF PACKAGING  
LEADERSHIP—AND THE REST IS HISTORY.**

Visit [comar.com](http://comar.com) or dial 1-866-655-7165 for more info.

## Oral Syringe Systems

Tested, approved, and used by major pharmaceutical companies, our oral syringes have a unique patented wiper design and meet quality standards that are unmatched in the industry.



### Accurate

User-friendly, frictionless operation  
Calibration lines vision-inspected for accuracy  
Pressure tested during assembly to ensure functionality

### Safe

One-piece plunger design eliminates the rubber grommet typically used, eliminating a choking risk  
Reduced extractables - latex free

### Durable & Reusable

Dishwasher safe  
Maintains print quality and seal integrity over extended usage

### Plunger Options

Material: HDPE  
Standard colors: red, yellow, blue, green, white

### Barrel Options

Material: polypropylene  
Stock colors: clear, amber tint  
Standard colors: white, red, yellow, blue, or green tint  
UV amber barrels meet USP 671 for light transmission

### Customization

Logos & lines - 360°, 1 or 2-color custom printing  
Colors - plungers, translucent barrels, tips caps, and stands  
Wrapping - bag/pouch, shrink wrap, wrapped with bottle adapter  
Tip options - for compatibility with other adapters and fittings

## Oral Syringes



Item Family	Description	Size (mL)	Tip
41-0129	1/2 mL Oral Syringe	.5	Short
41-0008	1 mL Oral Syringe	1	Short
41-0002	3 mL Oral Syringe	3	Long
41-0001	3 mL Oral Syringe	3	Short
41-0004	5 mL Oral Syringe	5	Long
41-0003	5 mL Oral Syringe	5	Short
41-0006	10 mL Oral Syringe	10	Long
41-0005	10 mL Oral Syringe	10	Short
41-0134	20 mL Oral Syringe	20	Short

## Tip Caps and Stands



Item Family	Description	Resin
14-0209	Oral Syringe Tip Cap	HDPE
14-0878	Oral Syringe Self-Righting Tip Cap	HDPE
14-0223	Oral Syringe No-Choke Tip Cap	HDPE


## Bottle Adapters



Item Family	Description	Style	Size (mm)	Resin
22-1230	DoseGuard®, Press-In Valved Bottle Adapter	Press-In	20	LDPE
22-0454	Press-In Bottle Adapter	Press-In	20	LDPE
22-0625	Vented Press-In Bottle Adapter	Press-In Vented	20	LDPE
22-0218	Press-In Bottle Adapter	Press-In	24	LDPE
22-0198	Press-In Bottle Adapter	Press-In	28	LDPE
22-0231	Press-In Bottle Adapter	Press-In	33	LDPE

### DoseGuard® Press-In Valved Bottle Adapter

Adds another level of protection for child-safe liquid medicine packaging and delivery.

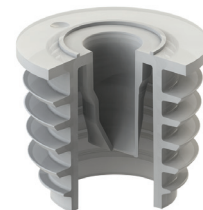

- Designed to function in industry standard bottles of HDPE, LDPE, PP and PET
- Adapter securely connects to the container finish and rests on top of the finish, this low profile allows for use with standard closures
- Features a patented valve design that remains closed until the tip of the oral syringe is inserted into the valve
- Valve opening is customizable to specific product viscosities and is easily adaptable to an open hole flow restrictor
- Single material one-piece construction reduces contamination risks

0279 OralSyringe 02.22

# 73 YEARS OF PACKAGING LEADERSHIP— AND THE REST IS HISTORY.

Visit [comar.com](http://comar.com) or dial 1-866-655-7165 for more info.